




# Intelligent Delivery...in a Flash



Latex Free • Fill Volume 2.0 cc  
Easy to Use: NO SALINE FILL REQUIRED

To learn more go to  
[www.activapatch.com](http://www.activapatch.com)


Distributed by:

Product of:  ActivaTek™



# Intelligent Delivery...in a Flash

**No Saline Needed**  
(Hydrogel ground electrode)

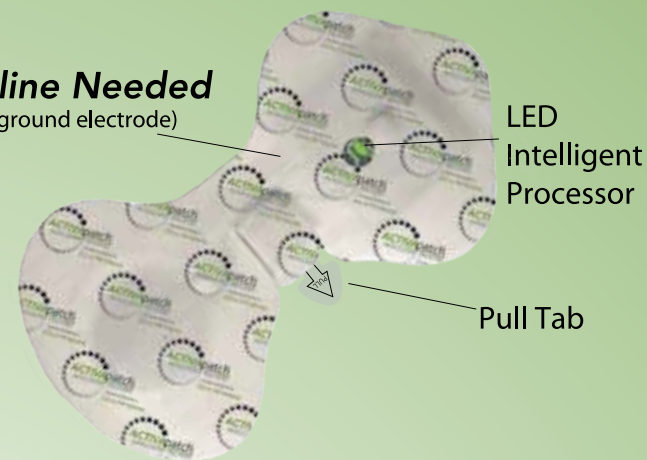

- . Intelligent
- . Effective
- . Fast
- . Accurate
- . Comfortable

See our interactive website  
[www.activapatch.com](http://www.activapatch.com)

Intelligent Delivery...in a Flash



# Smart Delivery System

## Effective

- Fully integrated SMART POWER system
- Constantly monitors skin resistance
- Detects problems with current continuity
- Indicates treatment mode
- No dose controller

## Convenient

- No external power source
- No saline fill
- Simply pull the tab to activate, fill the drug reservoir and apply directly to treatment site

## Fast & Accurate

- 80 mA\*min in 2.5 hrs  
40 mA\*min in 1.25 hrs
- Realtime power-off indicated by LED
- Most accurate delivery time; avg.  $\pm 2\%$  vs.  $\pm 37\%$  for other patches

## Comfortable

- Minimizes skin sensitivity
- Eliminates itching
- Thin and flexible
- No skin adhesive burns

Wireless Interactive Treatment	ActivaTek, Inc. ActivaPatch®
Power Source	6V
Interactive Intelligent Microprocessor	•
Interactive Treatment With Shut Off	•
Skin Resistance Indicator	•
Open Circuit Indicator	•
Memorizes Treatment Time	•
Dosage Options (40/80mA*min)	•
Treatment Time (80 mA*min)	2 hrs 30 min
80 mA*min Dosage	80mA*min $\pm 1.6$ mA*min avg.

### ActivaPatch® Specifications

Active Area	11.4 cm <sup>2</sup>
Dosage	40 or 80 mA*min
Skin Adhesive	Hypoallergenic
Maximum Buffer Capacity	Ag/AgCl Stabilizer
Fill Volume	2.0cc
Smart Power	6V with Intelligent Microprocessor
Packaging	6 ActivaPatch Kits, Alcohol Wipes, Patient Information Cards

Intelligent Delivery...in a Flash

A Product of: ActivaTek™

